

Brennen Ski

The non-carrier based instrument that makes meshing safe and easy.

The Brennen Skin Graft Mesher is a non-crushing, non-carrier based instrument that makes meshing safe and easy. This state-of-the-art mechanism meshes by “pinching” not cutting the graft, creating a uniform pattern and maintaining graft integrity.





Carrier-free

The mesher does not require plastic carriers, allowing you to mesh any graft length without being dictated by the length of a carrier. This also saves the time, trouble and cost of ordering disposable carriers.

Safe and easy to use


The Skin Graft Mesher is designed for convenience and ease of use. The smooth operation of the ratchet handle makes the meshing process effortless. The weight of the base provides added stability and holds the mesher firmly in place for unassisted operation.

The Skin Graft Mesher does not contain blades, eliminating the need for sharpening and reducing the concern for injury. The crossbar gives an extra element of steadiness when meshing and makes the unit easily transportable.

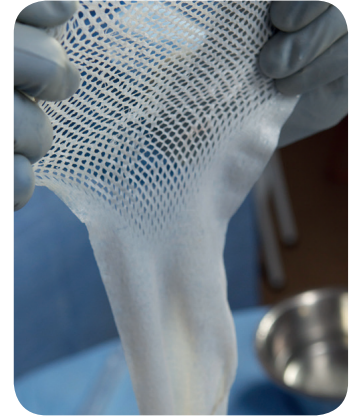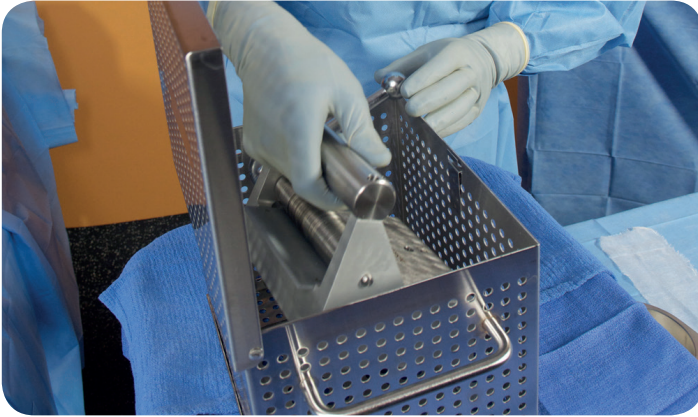

Skin graft mesher patterns


[representation using paper model]

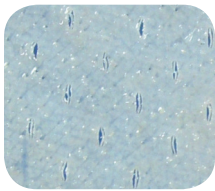


pie crust
(minimal perforation)




1:1

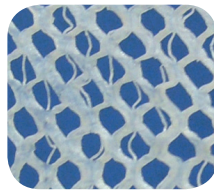


2:1

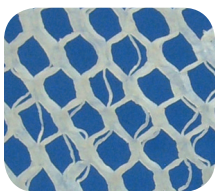

3:1

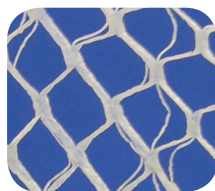

4:1

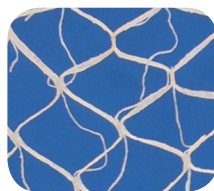

6:1

4.5 inch Skin Graft Mesher Assortment

Measures grafts up to 4.5 in (11.5 cm) wide

Art. No	Product Description
131500	Pie Crust (Minimal perforations)
131501	1:1 Expansion Ratio
131502	2:1 Expansion Ratio
131503	3:1 Expansion Ratio
131504	4:1 Expansion Ratio
131506	6:1 Expansion Ratio
131510	Mesher Storage/ Sterilizer Box

6.5 inch Skin Graft Mesher Assortment

Measures grafts up to 6.5 in (16.5 cm) wide

Art. No	Product Description
131600	Pie Crust (Minimal perforations)
131601	1:1 Expansion Ratio
131602	2:1 Expansion Ratio
131603	3:1 Expansion Ratio
131504	4:1 Expansion Ratio
131605	6:1 Expansion Ratio