

Tubifast™ 2-WAY STRETCH™

Lightweight, tubular
dressing retention


- **Holds dressings securely** without constriction or compression

- **Tubular design**
– no requirement for pins or tape.

- **Quick and easy to use:** cut to size and ease into place over the dressing for an even, non-constrictive fit

- Quickness and ease of use means the **risk of cross-contamination is reduced**

- **Less fray**, so there is less risk of fibres in the wound²

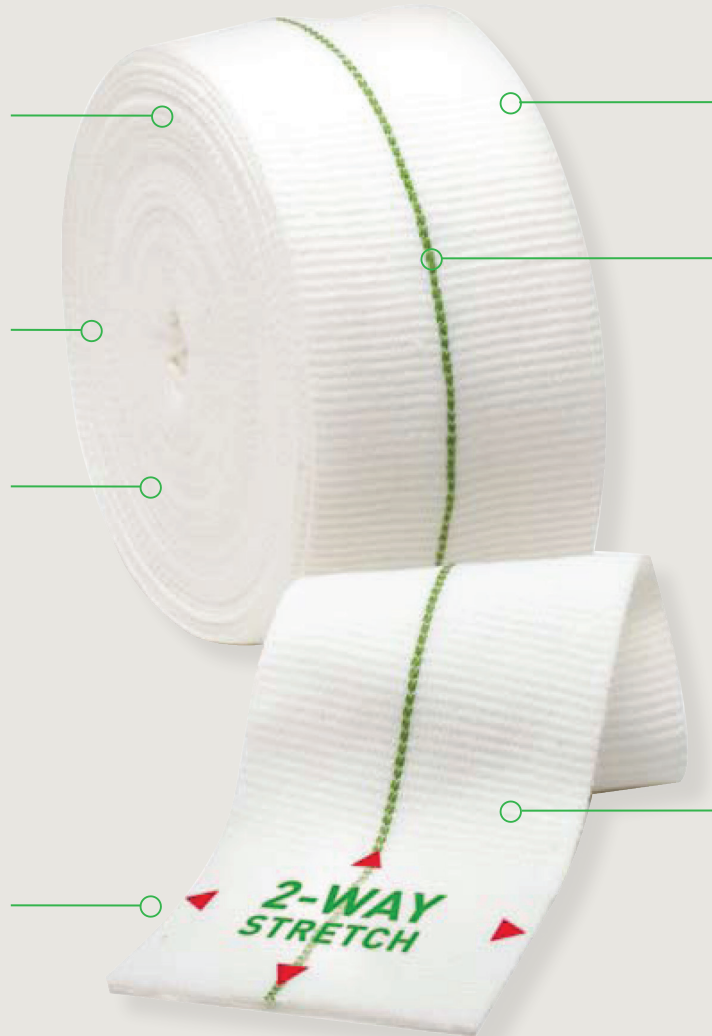

- **Cost effectiveness** with minimal waste: cut to size as required

- **Quick reference colour coding** to fit everything from small limbs to adult trunks


- Tubifast 2-Way Stretch's **light, radial and longitudinal elasticity** allows patients complete freedom of movement¹

- Tubifast™ 2-Way Stretch™ is the world's first and only lightweight tubular bandage with both radial and longitudinal stretch.¹ This unique mode of action allows patients complete freedom of movement, with added comfort. Tubifast™ 2-Way Stretch™ is ideal for dressing retention or as a skin covering for any part of the body. It can also be used for patch wrapping and as an under-cast stockinette. Tubifast™ 2-Way Stretch™ holds dressings securely in place without constriction or compression. It is quick and easy to use, with less risk of fraying² than traditional tubular bandages.


Tubifast™
2-WAY STRETCH™


Mölnlycke®


Directions for use

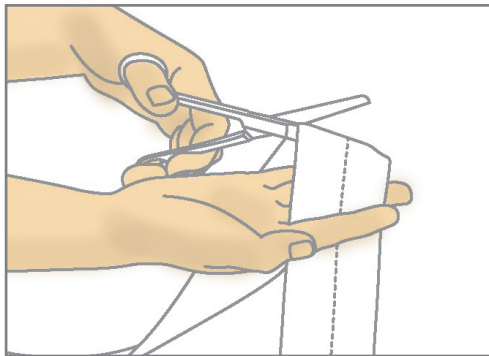

1. Choose the correct size, using the quick reference colour coding, and cut to length.

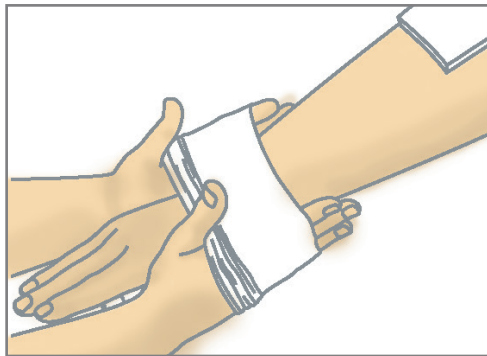

2. Stretch Tubifast 2-Way Stretch and ease gently into place over the affected area.

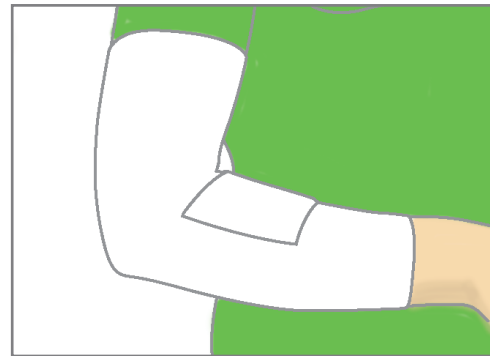

3. Position Tubifast 2-Way Stretch to cover the dressing.

How Tubifast 2-Way Stretch works

Tubifast 2-Way Stretch is manufactured from viscose with very fine elastane threads knitted into the fabric radially and longitudinally, to provide light elasticity.

Tubifast 2-Way Stretch holds dressings securely, without constriction or compression. No pins or tapes are necessary. Its two-way stretch allows patients complete freedom of movement and added comfort.

Tubifast 2-Way Stretch is available in a range of sizes designed to use on any part of the body. Each size is identified by a coloured line knitted into the bandage. It's quick and easy to use: cut to size and stretch over the dressing for an even, non-constrictive fit.

Benefits of Tubifast™

- Tubifast 2-Way Stretch holds dressings securely, without constriction or compression.
- Its radial and longitudinal stretch allows patients complete freedom of movement.
- A variety of lengths are available, ensuring cost-effectiveness by minimizing waste.
- Quick and easy to use: simply cut to size and ease gently into place over the dressing for an even, non-constrictive fit.
- Available in a range of quick reference, colour-coded sizes to fit everything from small limbs to adult trunks.
- It has less fray so there is less risk of fibres in the wound².

Fibre content

Viscose, Elastane and Polyamide

Washing Instructions



Notice: For Mölnlycke licensed product details including indication and precaution, please refer to www.molnlycke.ca

Natural rubber latex is not a constituent of Tubifast 2-Way Stretch or its packaging.

¹ Published data shows that Tubifast 2-Way Stretch exceeds the minimum requirement of longitudinal stretch for the actions of sitting, bending and flexing. Physical testing of textiles, B. P. Saville, Woodhead Publishing, 1999.

² Data on file.

Mölnlycke
1-800-494-5134
www.molnlycke.ca

Mölnlycke Health Care AB, Box 13080, SE-40252 Göteborg, Sweden. Phone +46317223000.
The Mölnlycke Health Care, Tubifast™ and 2-Way Stretch™ names and respective logos are registered globally to one or more of the Mölnlycke Health Care Group of Companies.
© 2019 Mölnlycke Health Care AB. All rights reserved.

Indications*

For dressing retention, wet and dry patch wrapping. Select the correct size and cut to length, open the bandage and stretch over the affected area. A double layer can be applied for extra security.

Precautions*

- If soiled with bodily fluids, discard appropriately and replace. Avoid naked flames and ignition sources, especially when used in conjunction with paraffin based products. If condition worsens discontinue use and seek medical advice.

Tubifast 2-Way Stretch Assortment

Art. No.	Product	Width x Length
RED LINE 2434	Small limbs	3.5cm x 10m
GREEN LINE 2436	Small and medium limbs	5cm x 10m
BLUE LINE 2438	Large limbs	7.5cm x 10m
YELLOW LINE 2440	Extra-large limbs, heads, children's trunks	10.75cm x 10m
PURPLE LINE 2444	Adult trunks	20cm x 10m

Tubifast 2-Way Stretch is packaged in single packs




Mölnlycke®
CEWC01521905